

# Bright medico

(Never stop learning)

**BRIGHT**  
**MEDICO**

# **LAB MANUAL PHARMACY**

**1ST YEAR**

**PHARMACOGNOSY**



[www.Brightmedico.com](http://www.Brightmedico.com)

@brightmedico

WEBSITE - [WWW.BRIGHTMEDICO.COM](http://WWW.BRIGHTMEDICO.COM)

Get in touch with us - search (BRIGHT MEDICO)

**Bright medico**  
(Never stop learning)

## **PRACTICAL LAB MANUAL**

### **PHARMACOGNOSY**

**D. Pharm I<sup>st</sup> Year**

# Bright medico

(Never stop learning)

Sr. No	Name of Experiment
1.	To study the compound microscope
2.	To study the macroscopic and microscopic study of the Senna Leaf.
3.	To study the macroscopic and microscopic study of the Datura Leaf.
4.	To study the macroscopic and microscopic study of the ISAPGOL
5.	To study the macroscopic and microscopic study of the CARDAMOM
6.	To study the macroscopic and microscopic study of the FENNEL
7.	To study the macroscopic and microscopic study of the DILL
8.	To study the macroscopic and microscopic study of the CORIANDER
9.	To study the macroscopic and microscopic study of the Clove
10.	To study the macroscopic and microscopic study of the CINNAMOM
11.	To study the macroscopic and microscopic study of the GINGER
12.	To study the test of the AGAR:
13.	To study the test of the acacia
14.	To study the test of the TRAGACANTH
15.	To study the test of the GELATIN
16.	To study the test of the STARCH
17.	To study the test of the HONEY
18.	To study the test of the CASTOR OIL


# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -1

**AIM:** To study the compound microscope

**REQUIREMENTS:** Compound Microscope, Permanent slide (or any other specimen).

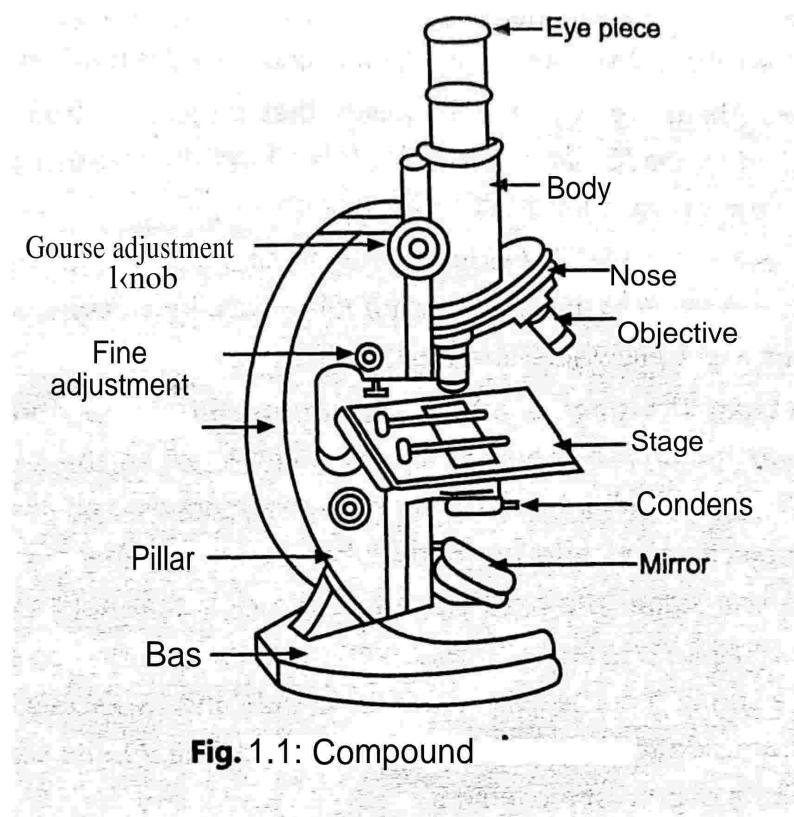

**Fig. 1.1: Compound**

### Theory:

The microscope is one of the most commonly used instruments in medical, paramedical and clinical laboratories. It is used to study Cell Morphology, Histology, Histopathology and Microbiology. A microscope helps us to see microscopic objects that are too small and invisible to the naked eye.

### Description of Compound Microscope:

The Compound microscope has the following main parts

1. The supporting system.

# Bright medico

## (Never stop learning)

2. The focusing system.
  3. The optical or magnifying system.
  4. The illumination system including:
    - a) Source of light
    - b) Mirror
    - c) Condenser
1. **The support system:** It is a framework to which various functional units are attached. It consists of the following:
- a) **Base:** it is a heavy metallic, 'U' shaped or 'horseshoe' shaped base with supports the microscope on work table and provides maximum stability.
  - b) **Pillars:** There are two upright pillar that project up from the base and are attached to the 'C' shaped handle. This allows the microscope to be tilted at a suitable angle for comfortable observation.
  - c) **Body tube:** It is 16-17 cm long cylindrical tube fitted at the upper end of handle which is vertical or at an angle through which light passes via the eye piece to the observer's eye visualizing the image
  - d) **The stage:** The stage is a square platform with an aperture in its centre and fitted to the limb below the objective lenses. When the slide is placed on it, converging rays of light emerging from the condenser passes through slide and then objective lens into the body tube. It can be either the fixed stage or the mechanical stage. The fixed stage has two clips that hold the slide in position. The mechanical stage has a calibrated metal frame fitted on right side of stage. It has a spring mounted clip to hold the slide and two screw heads to move the slide from side to side, forward and backward. The Vernier scale is also attached to indicate degree of movement

**The focusing system:** The focusing system consists of coarse and fine adjustment and screw heads are used for raising and lowering body tube for proper focusing the slide. The coarse adjustment moves the focusing system up or down through a large distance via a rack. The fine adjustment works in same way which requires several rotations to move the tube through a small distance. It is employed for accurate focusing.

**The optical or magnifying system:** It consists of body tube, eyepiece and the nosepiece. The body tube is the present between upper end of objective and eyepiece. The eyepiece fits into top of body tube. They can be 5X, 6X, 8X, 10X or 15X. Each eyepiece has two lenses; the eye lens at the top and field lens at the bottom. The field lens collects divergent rays, passes through eye lens to further magnify the image. The Nosepiece has two parts; the fixed nosepiece and the revolving nosepiece. The fixed nosepiece holds the revolving nosepiece that carries interchangeable objective lenses. The objective lenses are spring loaded objectives of different magnifying powers.

Different types of the objective lenses are low power objective or 10X, high power objective or 45X, oil immersion objectives Of 100X and scanning objective 3X.

**Magnification:** Magnification is the ability to make small objects seen larger, such as making a microscopic organism visible. The objective lenses magnify the images as stated below:

Low power objective  $10X = 10 \times 10 = 100$  times. High power objective  $45X = 45 \times 10 = 450$  times.

oil immersion objective (100X)  $= 100 \times 10 = 1000$  times.

oil immersion (100X) objective has very small aperture and deep focusing position i.e. 1 mm from the

# Bright medico

## (Never stop learning)

slide. The light rays coming from the slide (denser medium) are refracted by thin layer of air (rarer medium) away from the small aperture of the objective and results in faint image. If some other medium like cedar wood oil, paraffin or glycerin having same refractive index as that of glass is added on the slide, it removes the thin layer of air and forms a continuous medium. This avoids the refraction of light rays and results in sharp image.

**2. The illumination system:** The microscope will function only when proper illumination or lightning is provided. The illumination system is to provide uniform, soft and bright illumination.

The illumination system consists of:

(a) **Source of light:** It may be external (natural day light, electric lamp or tube light) or internal (electric in-built light source).

(b) **A condenser:** It is a system of lenses filled as short cylinder mounted below the stage.

(r) **Mirror:** A double sided mirror with one flat side and other concave is located below condenser and can be rotated in all directions. It focuses light rays into a solid cone of light onto the material under study and helps in resolving image.

(d) **Iris diaphragm:** Iris diaphragm is a thin opaque membranous structure fitted within condenser with a small lever on the side. The lever can adjust size of aperture of diaphragm and allows less or more light falling on slide,

### Procedure:

1. Examine the permanent slide/blood film/specimen first with naked eye.
2. Place the microscope on working table in an upright position, and raise the body tube approximately 7-8 cm above the stage. Put the slide on the stage and using the mechanical stage, bring the specimen over the central aperture.
3. Select the low magnification objective (10X). ‘
4. Select and adjust the mirror (plane or concave) so that the light shines on the specimen 5. Adjust the condenser well down, and partly close the diaphragm to cut down excess light.
6. Looking from the side, and using the coarse adjustment, brings the body tube down so that the low power lens is about 1 cm above the slide. Look into eyepiece and gently raise the tube till the slide comes into focus.
7. Then choose the area of interest for viewing it under higher magnifications.
8. For focusing under high magnification, simply rotate the nosepiece so that the high magnification objective (45X) ‘clicks’ into position. Raise the condenser to mid position and open the diaphragm to admit enough light. Use fine adjustment as required.
9. For focusing under oil immersion objective (100X), raise the body tube 8-10 cm above the slide. Place a drop of cedar wood oil, paraffin or glycerin on & the slide. Looking from the side bring down the objective till it just enters the oil drop. Use other adjustment as required.

# Bright medico

(Never stop learning)


## EXPERIMENT NUMBER -2

**AIM:** To study the macroscopic and microscopic study of the Senna Leaf.

**SYNONYMS:** Senai ki patti, Sonamukhi, Tinnevelly senna.


**BIOLOGICAL SOURCE:** It consists of dried leaflet of cassia angustifolia Vahl. Belonging to family Leguminosae and contained not less than 2.0% of glycosides calculated as Sennoside B.

**Macroscopy:**

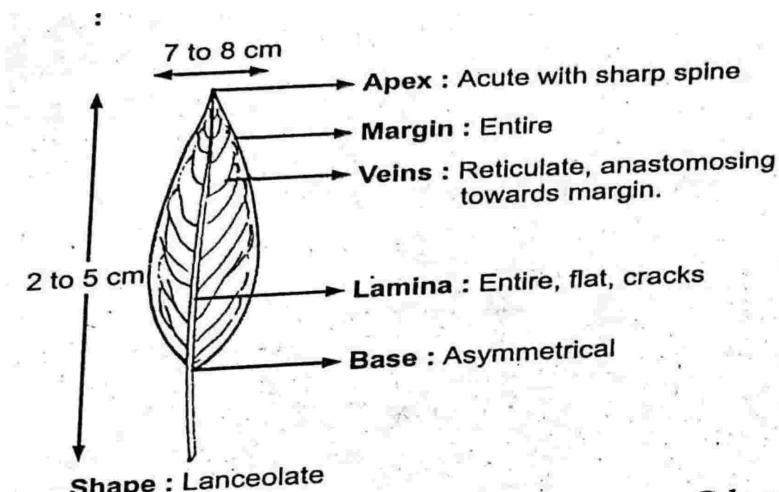

**Organoleptic CHARACTERS:** Colour: Yellowish, Odour: slight, Taste: Mucilaginous, slightly bitter.

**EXTRA FEATURES:** Surface: Isobilateral, thin, pubescent (hairy) with trichomes on both surfaces.

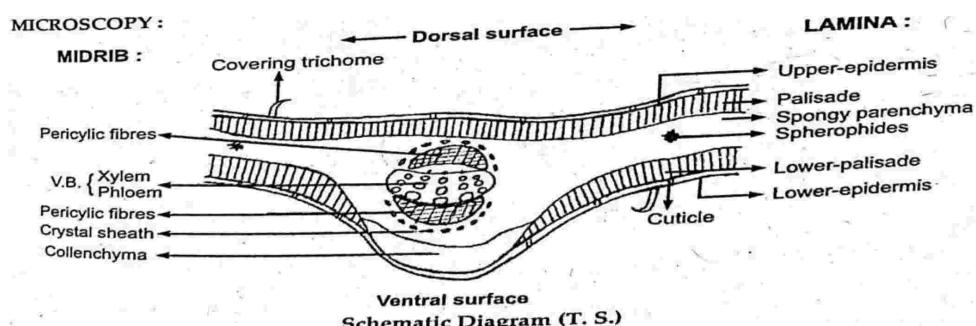


# Bright medico

(Never stop learning)


## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:

SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ HCl(1:1)	Conc.	RED/PINK Lignified tissues: xylem(vascular bundle), sclerenchyma (pericycle)
2	Ruthenium red		RED/PINK Mucilaginous cells of epidermis
3.	Sudan red III	pink	Cutin/cuticle
4.	Acetic acid	Crystal insoluble	Calcium oxalate
5.	Hydrochloric acid	Crystal soluble	Calcium oxalate
6	Sulphuric acid (60%)	Crystal soluble	Calcium sulphate needles

## SURFACE PREPARATION : (For procedure see the topic LEAVES)

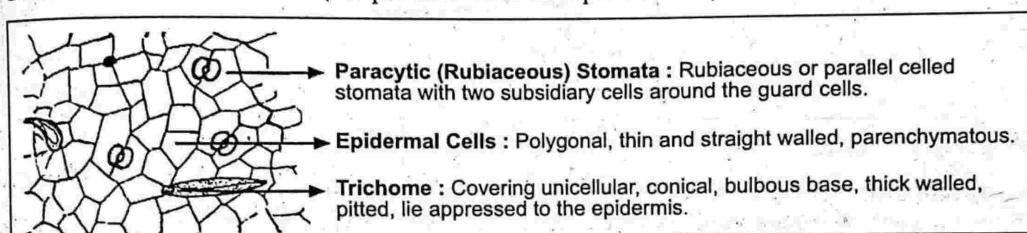


## MICROSCOPIC CHARACTERISTICS OF THE POWDER DRUGS:

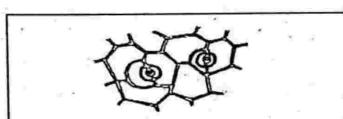

**Epidermal Cells :** Polygonal, straight walled, epidermal cells, with paracytic stomata (Rubiaceous)


**Paracytic Stomata :** See the description above.

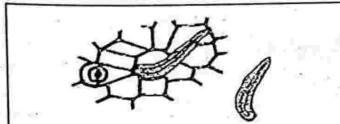


**Covering Trichomes :** Unicellular, thick warty walls, acute apex, bulbous base, narrow lumen, conical shape.

Length = 70 → 260  $\mu$

Width = 12 → 18 → 25  $\mu$

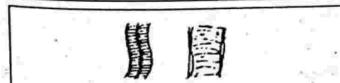

**Xylem Vessels :** Angular thickening, lignified.

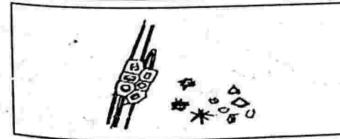

**Calcium Oxalate :** Crystals isolated or in parenchymatous cells. Very abundant, occur as prisms and also as cluster crystals.

# Bright medico

(Never stop learning)

## **CHEMICAL TEST:**

**BRONTRAGER TEST FOR ANTHRAQUINONE:** boil the leaves with dil.sulphuric acid. Filter and cool the filtrate. Add immiscible organic solvent layer in another test tube. add strong ammonia solution, shake slightly and keep the test tube aside , lower ammonical layer shows pink/ red colour.

**Chemical constituents:** Mainly anthraquinone glycosides: Sennosides A, B, C, D.

**Uses:** Irritant purgative.

# Bright medico


(Never stop learning)

## EXPERIMENT NUMBER -3

**AIM:** To study the macroscopic and microscopic study of the Datura Leaf.


**Synonym:** Datura herb

**Biological Source:** Datura mainly consists of the dried leaves and flowering tops of the Datura metel Family : Solanaceae. It contains not less than 0.2% of total alkaloids calculated as Hyoscyamine.

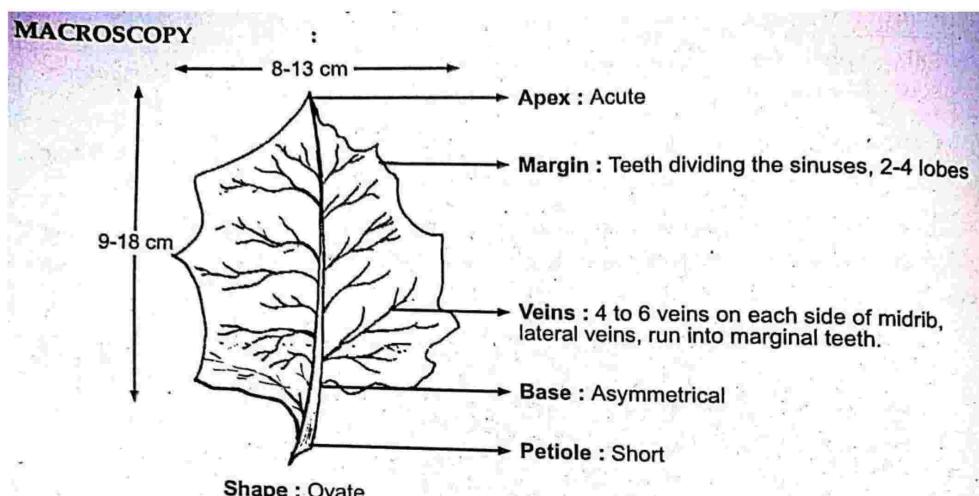

**ORGANOLEPTIC CHARACTERISTIC:** Colour: Pale green, odour: disagreeable characteristic Taste: unpleasant bitter.

**EXTRA FEATURES:** Texture: Thin and minutely hairy upper epidermis darker than lower, midrib prominent on lower surface.

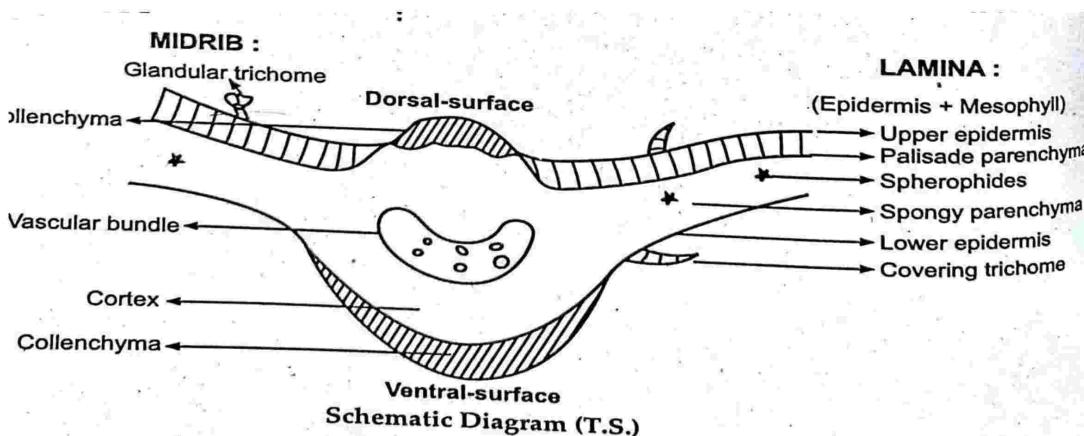


# Bright medico

(Never stop learning)

## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:


SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ HCl(1:1)	Conc.	RED/PINK Lignified tissues: xylem(vascular bundle))
2.	Dill Acetic acid	insoluble	Calcium oxalate
3.	Sulphuric acid (60%)	Crystal soluble	Calcium sulphate needles

### SURFACE PREPARATION : (For procedure see the topic LEAVES)

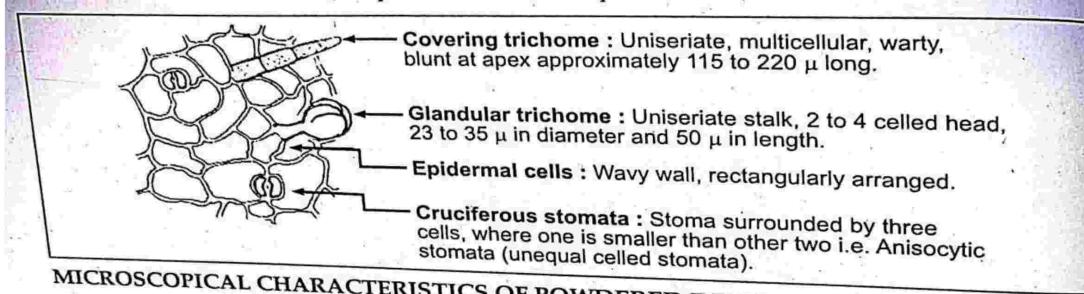

### MICROSCOPICAL CHARACTERISTICS OF POWDERED DRUG :


#### STOMATA AND EPIDERMAL CELLS

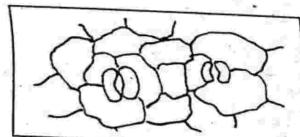

**Cruciferous Stomata (Anisocytic) :** See above

**Epidermal cells :** Slightly straight in upper epidermis, wavy walled in lower epidermis.

#### TRICHOMES

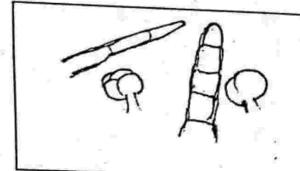


##### Covering Trichome :

Uniseriate, multicellular, wavy covering trichomes with blunt apex.  
Length= approximately 115 to 200  $\mu$ .

##### Glandular Trichome :


One celled stalk and 2-4 celled head,  
23 to 35  $\mu$  in diameter and 50  $\mu$  in length.

#### LAMINA / MESOPHYLL

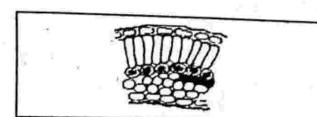

**Palisade Parenchyma :** Elongated, compactly arranged, Cal. oxalate crystals (spherophides), vascular element, Spongy parenchyma.

#### CALCIUM OXALATE CRYSTALS

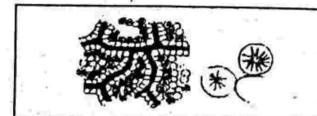

**Spherophide Crystals :** Occur in spongy parenchymatous cells and scattered in the powder.

# Bright medico

(Never stop learning)

## **CHEMICAL TESTS:**

1. Vitali Morin reaction: The tropane alkaloid is treated with fuming nitric acid, followed by evaporation to dryness. Addition of methanolic potassium hydroxide solution to an acetone solution of nitrated residue, Violet colour is developed.
2. On addition of silver nitrate solution to solution of hyoscine hydrobromide, yellowish white precipitate is formed which is insoluble on nitric acid, but soluble in dil. Ammonia solution.

**CHEMICAL CONSTITUENTS:** Up to 0.5% of total alkaloids Scopolamine (hyoscine) is the main alkaloid. Hyoscyamine and atropine are present in minor quantities.

**USES:** Parasympatholytic with anticholinergic and CNS depressant effect. Also used as mydriatic , antispasmodic and cerebral sedative. Hyoscine hydrobromide is used in motion sickness, gastric or duodenal ulcer.

# Bright medico


(Never stop learning)

## EXPERIMENT NUMBER -4

**AIM:** To study the macroscopic and microscopic study of the ISAPGOL


**Synonym:** Isapghula, Isapgol, Indian Psyllium

**Biological Source:** Dried seeds of *Planago ovata* Forsk Family: Plantaginaceae.

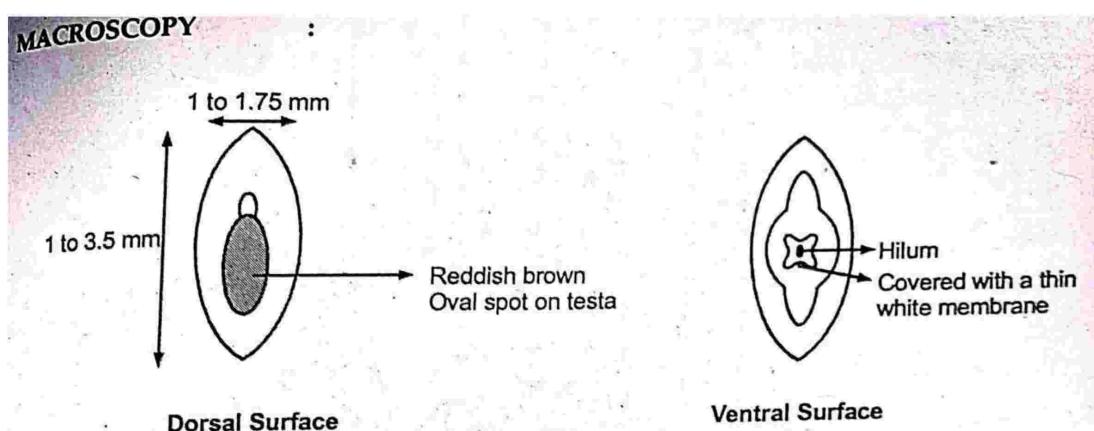

**ORGANOLEPTIC CHARACTERS:** Colour: Pinkish grey to brown, odour : none, taste : Mucilaginous.

**EXTRA FEATURES:** Testa is hard, translucent and smooth, Weight of 100 seeds = 0.15-0.19 gms, Swelling factor= 10.25-13.50

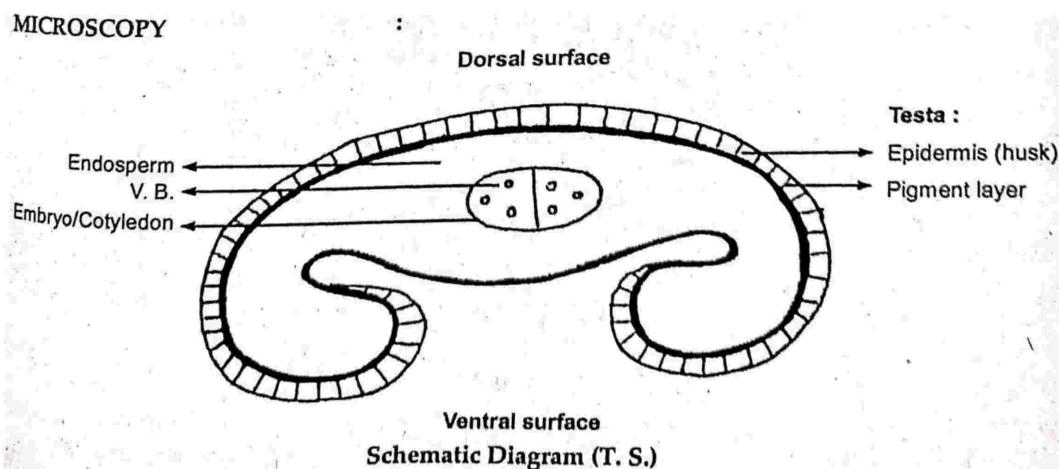


# Bright medico

(Never stop learning)


## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:

SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ HCl(1:1)	Conc.	RED/PINK Lignified tissues: xylem(vascular bundle))
2.	Ruthenium Red	Red	Mucilage present in epidermis
3.	Alcoholic Picric acid	Yellow	Aleurone grains present in the cells of endosperm and embryo
4	Sudan Red III	Red	Oil globules present in the cells of endosperm and embryo

### MICROSCOPICAL CHARACTERISTICS OF POWDERED DRUG :

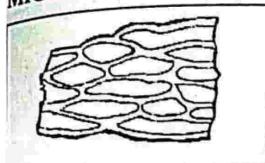


**Epidermis :** Thick walled transparent, polyhedral to tangentially elongated cells, containing mucilage. The very occasional starch granules, which are present in some of the epidermal cells and may be found embedded in the mucilage. They are very small and simple or compound.

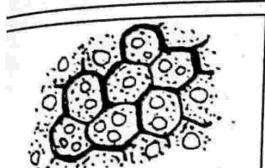

**Endosperm :** Thick walled cells with numerous pits. Oil globules show red colour with Sudan Red III.

Aleurone grains show yellow colour with alcoholic picric acid solution.

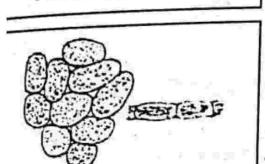

**Pigment Layer :** Yellow pigment layer.

### CHEMICAL TESTS:

TEST	OBSERVATION	INFERENCE
Seeds on slide + water	Zone of mucilage is formed around each seed	Mucilage present
Powder+ Water	Forms jelly like mass	Mucilage present

**CHEMICAL CONSTITUENTS:** Mucilage 10%, carbohydrates, fixed oil, proteins.

**USES:** demulcent, laxative in chronic constipation. Mucilage is used as a coating materials binding agent etc.

# Bright medico


(Never stop learning)

## EXPERIMENT NUMBER -5

**AIM:** To study the macroscopic and microscopic study of the CARDAMOM


**Synonym:** Chotti - Ilaychi

**Biological Source:** It consists of dried, nearly ripe fruit of *Elettaria cardamomum* Var. Family: Zingiberaceae. Seeds contain of less than 4% of volatile oil.

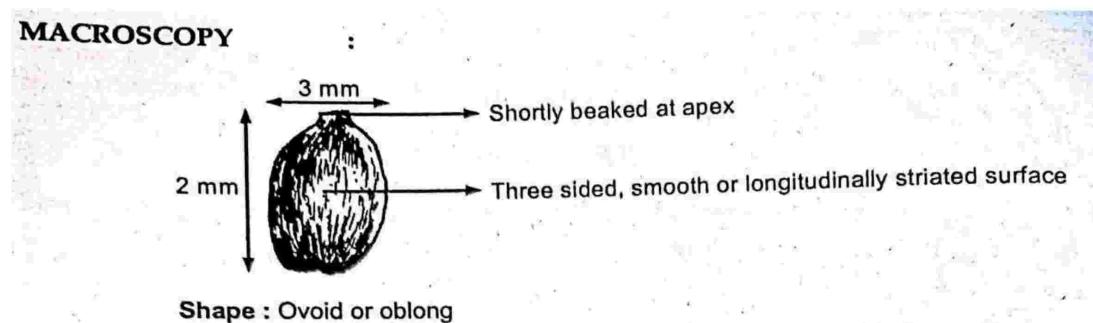

**ORGANOLEPTIC CHARACTERS:** Colour: Pericarp.: Green to pale buff, Seeds: pale to reddish brown , odour: strongly aromatic and characteristic, taste: strongly aromatic.

**EXTRA FEATURES:** fruits contain three chambers, each consists of two rows of seeds. Seeds are enclosed in membranous arillus. Each chamber contains six to ten seeds in two rows.

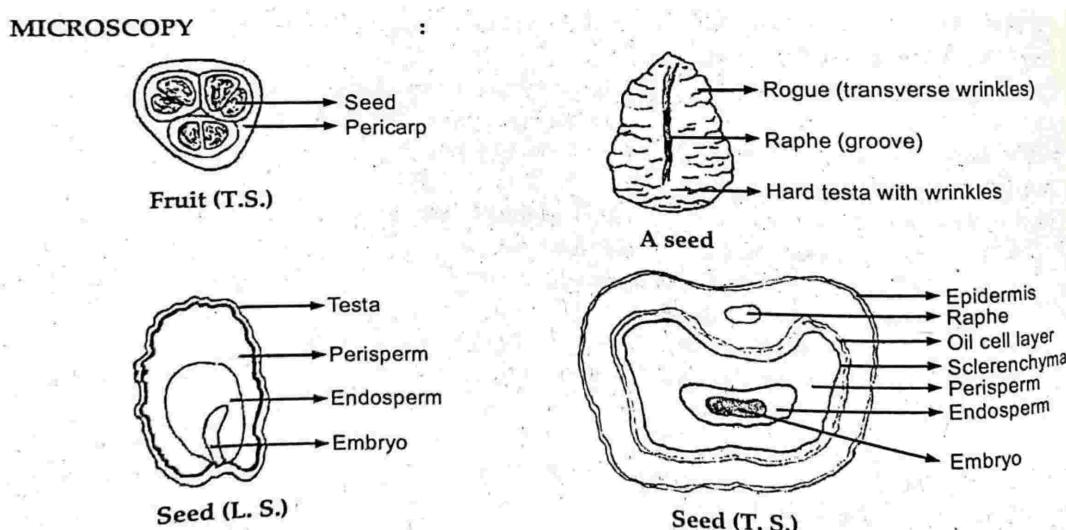


# Bright medico

(Never stop learning)

## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:

SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ Conc. HCl(1:1)	RED/PINK	Lignified tissues: xylem(vascular bundle))
2.	SULPHURIC ACID (60%)	RED	Calcium oxalate
3.	Dil. Iodine solution	blue	Starch present in perisperm
4	Sudan Red III	red	Oil globule
5	Alcoholic picric acid	yellow	Aleurone grains


## MICROSCOPICAL CHARACTERISTICS OF POWDERED DRUG :

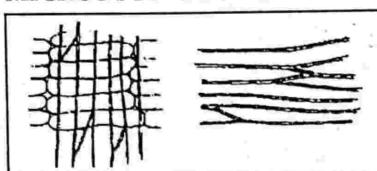

**Epidermis :** Epidermal cells are straight walled, colourless to light yellow, parallel to each other, may be associated with oil cells.



**Stone cells :** Lignified, red to orange coloured fragments, cylindrical and elongated in side view and polygonal in surface view.

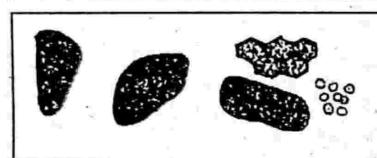

**Perisperm :** Thin walled cells with calcium oxalate crystals and starch grains.

**Starch grains :** 2 to 4 to 6  $\mu$  in diameter.

**Calcium oxalate :** 10 to 12  $\mu$  long. Scattered as well as in the cells of perisperm.

**CHEMICAL CONSTITUENTS:** **Volatile oils:** Cineol, limonene, borneol, L-terpenol, starch, fixed oil, silica.

**USES:** carminative, stimulant, aromatic, flavouring agent


# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -6


**AIM:** To study the macroscopic and microscopic study of the  
**FENNEL** Synonym: Bari sauf

**Biological Source:** Consists of dried ripe fruit of cultivated species *Foeniculum vulgare*. It contains not less than 1.4% of volatile oil. **Family:** Umbelliferae.

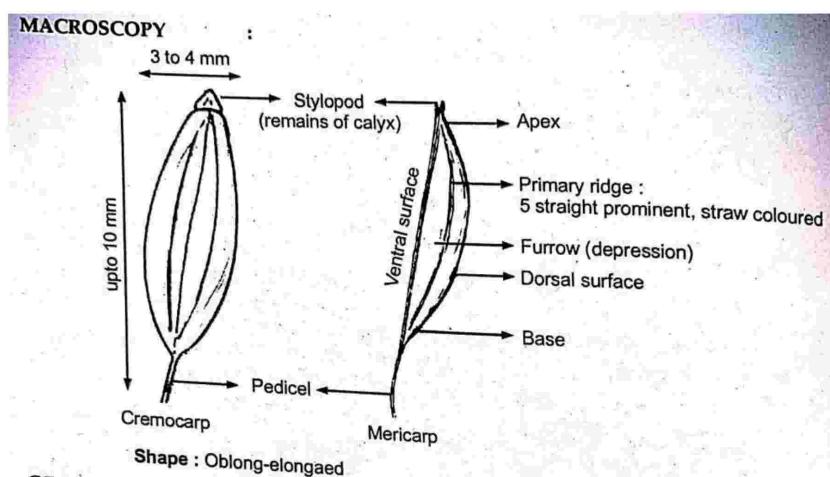

**ORGANOLEPTIC CHARACTERS:** Colour: Greenish or yellowish brown **Odour:** sweet aromatic and characteristic, **Taste:** Sweet, mucilaginous, agreeable, aromatic and characteristic.

**EXTRA FEATURES:** fennel exhibits **cremocarp**. It consists of two equal portions called as **Mericarps**, connected by central stalks called as **carpophores**.

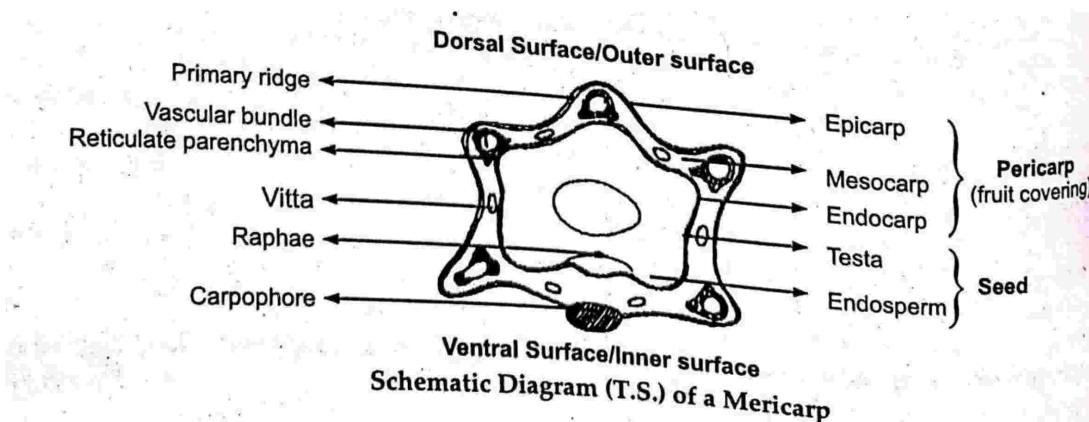


# Bright medico

(Never stop learning)

## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:

SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ HCl(1:1)	Conc.	RED/PINK Lignified tissues: xylem(vascular bundle))
2	Sudan Red III	red	Oil globule in the cells of endosperm and cuticle
3	Alcoholic picric acid	yellow	Aleurone grains

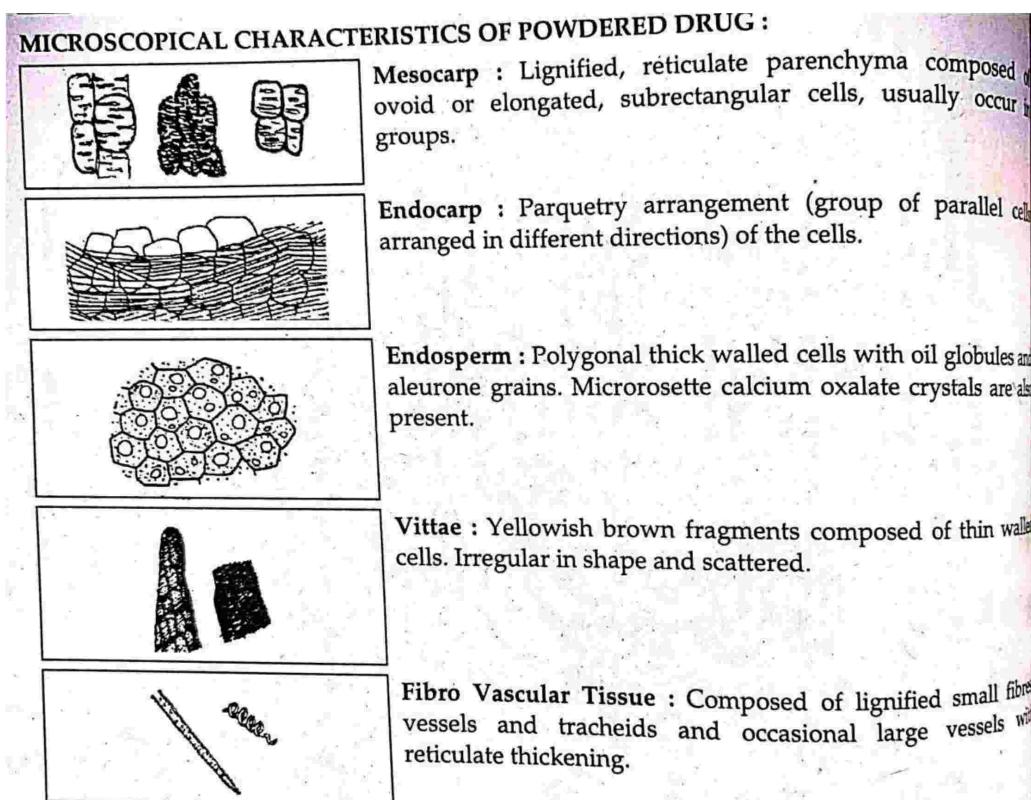

**CHEMICAL CONSTITUENTS:** **Volatile oils:** (4-6%), **Anethole** 50-60% of vol. oil, **d-fenchone** 10% of volatile oil, Fixed oil 12-18%, Proteins (14-22%)

**Uses:** carminative, respiratory stimulant, aromatic.

# Bright medico


(Never stop learning)

## EXPERIMENT NUMBER -7

**AIM:** To study the macroscopic and microscopic study of the DILL


**Synonym:** *Fructus anethi*, European dill, *Anethum*

**Biological Source:** Consists of dried ripe fruit of *Anethum graveolens*. It contains not less than 2.5% of volatile oil. **Family:** Umbelliferae.

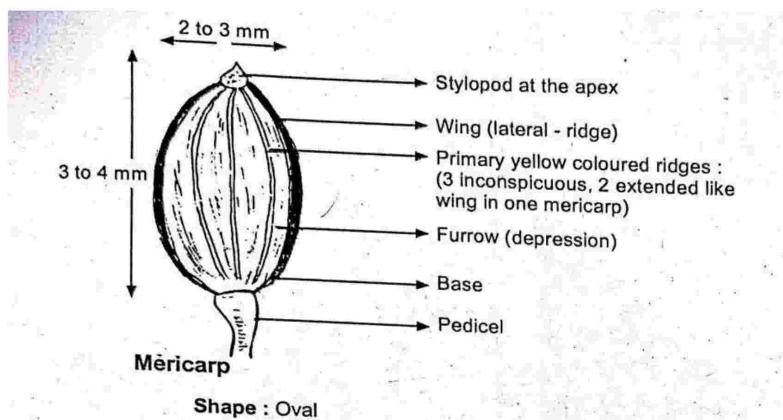

**ORGANOLEPTIC CHARACTERS:** Colour: chocolate brown, Odour: Aromatic, spicy and characteristic, Taste: Spicy, aromatic, and characteristic.

**EXTRA FEATURES:** normally separated mericarps, dorsally compressed. Orthospermous fruits.

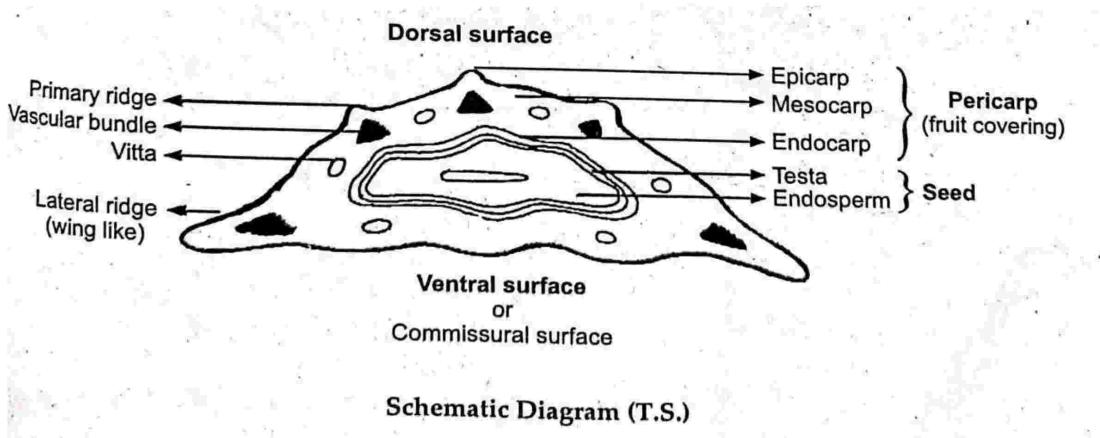


# Bright medico

(Never stop learning)

## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:

SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ HCl(1:1)	Conc.	RED/PINK Lignified tissues: xylem(vascular bundle))
2	Sudan Red III	red	Oil globule in the cells of endosperm and cuticle
3	Alcoholic picric acid	yellow	Aleurone grains

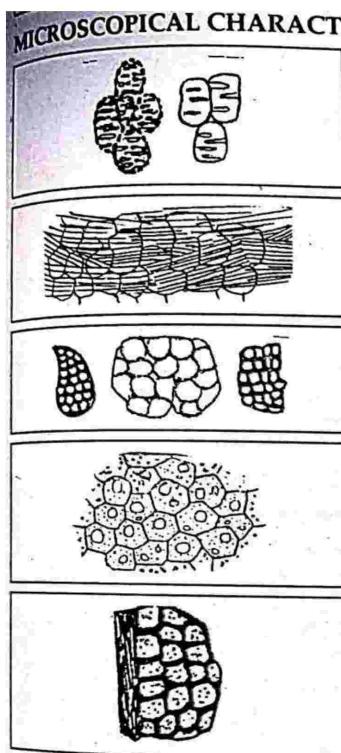

**Mesocarp :** Lignified parenchymatous cells of mesocarp showing reticulate thickenings.  
Groups of these cells are usually found associated with the fibrovascular tissue.

**Endocarp :** Parquetry arrangement is seen in surface view.

**Vittae :** Yellowish fragment with fine cracks, looks like broken glass. Composed of thin walled cells.

**Endosperm :** Thick walled polygonal cells, with oil globules, aleurone grains and micro rosette crystals of calcium oxalate.

**Sclerides :** Stone cells with rectangular, pitted walls. Lignified.

**CHEMICAL CONSTITUENTS:** volatile oils 3-4%, 50-60% carvone , 5% dihydrocarvone , D-limonene and phellandrene etc, fixed oil protein.

**USES:** Aromatic, stimulant, Antiseptic.

# Bright medico


(Never stop learning)

## EXPERIMENT NUMBER -8

**AIM:** To study the macroscopic and microscopic study of the CORIANDER


**Synonym:** Dhania (hindi)

**Biological Source:** Consists of dried ripe fruit of *Coriandrum sativum*. It contains not less than 0.3% of volatile oil. **Famil :** Umbelliferae.

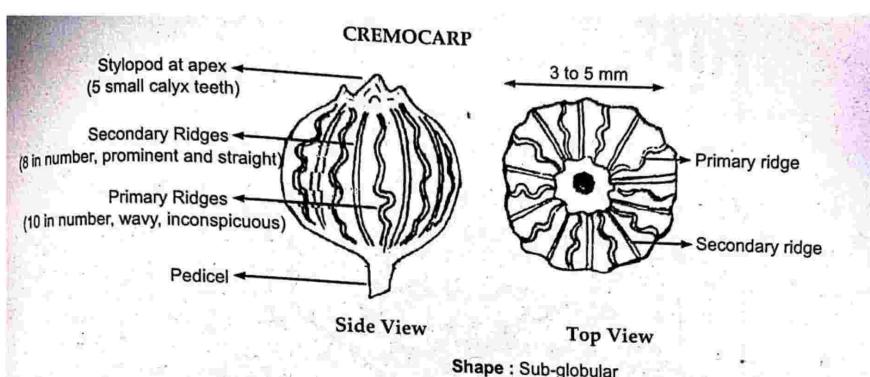

**ORGANOLEPTIC CHARACTERS:** Colour: Brownish yellow, Odour: Aromatic, taste: Spicy and characteristic.

**EXTRA FEATURES:** Cremocarps, caelospermous fruit.

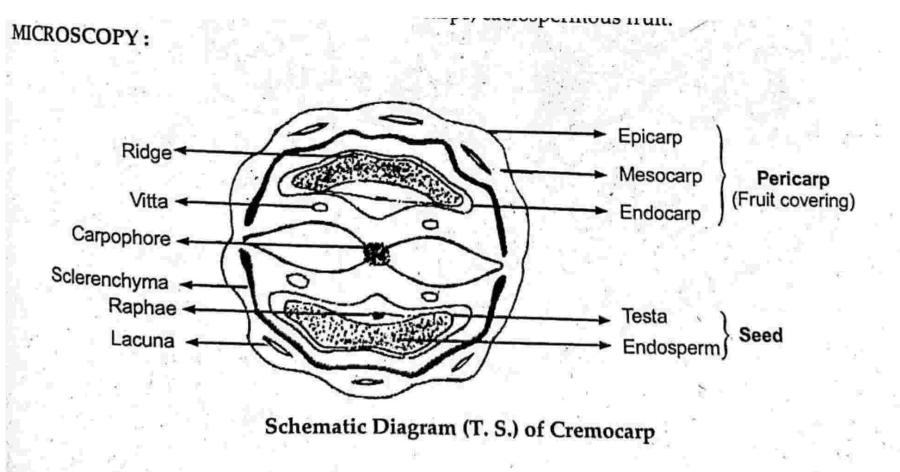


# Bright medico

(Never stop learning)


## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:

SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ HCl(1:1)	Conc.	RED/PINK Lignified tissues: xylem(vascular bundle))
2	Sudan Red III	red	Oil globule in the cells of endosperm and cuticle
3	Alcoholic picric acid	yellow	Aleurone grains


## MICROSCOPICAL CHARACTERISTICS OF POWDERED DRUG:

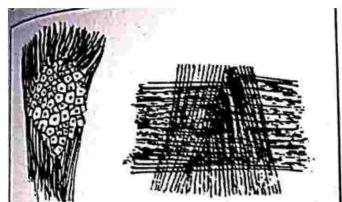


**Sclerenchymatous layer :** Groups of fusiform fibres running wavy; crossing each other, lignified.

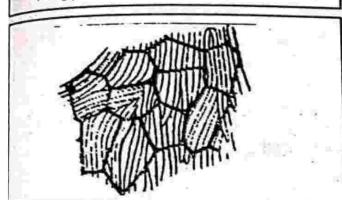

**Endocarp :** Parquetry arrangement of thin walled lignified cells with polygonal cells of mesocarp.

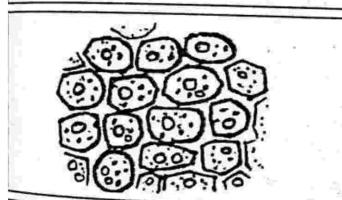

**Endosperm :** Polygonal parenchyma with aleurone grain and oil globules.  
Microrosette crystals of calcium oxalate in the cells.

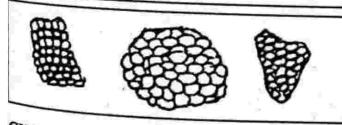

**Vittae :** Few yellowish brown fragments of vittae.  
Thin walled, small parenchymatous cells in group.

**Trichomes and lignified reticulate parenchyma are absent.**

**CHEMICAL CONSTITUENTS:** Volatile oil 0.2-1%: coriandrol (D-Linalool 60-70%).  
Terpene 20%, Fixed oil 13-20%, Proteins 17%.

**USES:** carminative, aromatic, stimulant, spice, flavouring agent.

# Bright medico


(Never stop learning)

## EXPERIMENT NUMBER -9


**AIM:** To study the macroscopic and microscopic study of the

**Clove Synonym:** Lavang (Hindi)

**Biological Source:** It consists of dried flower buds of *Eugenia caryophyllum*. It contains not less than 16% of Clove oil. **Family:** Myrtaceae.

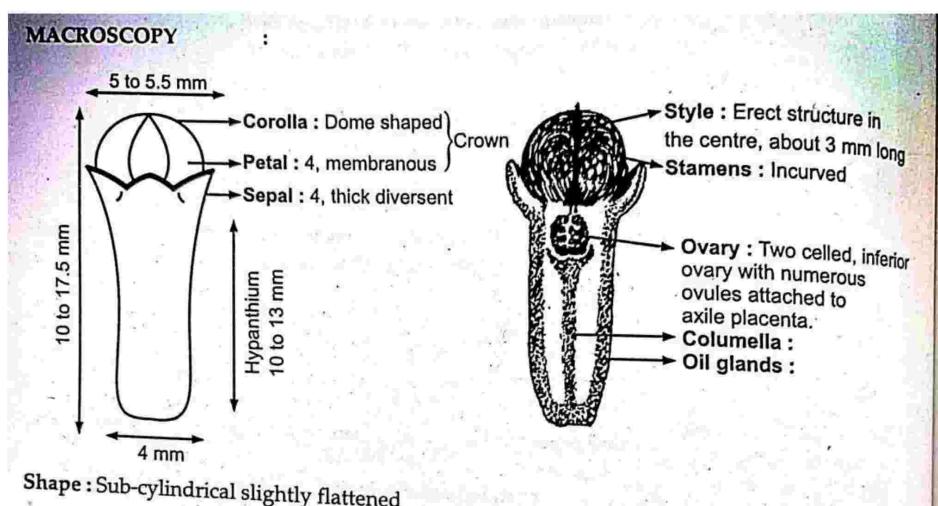

**ORGANOLEPTIC CHARACTERS:** Colour: Dark brown or crimson red, Odour: Aromatic, Taste; Spicy pungent followed by numbness.

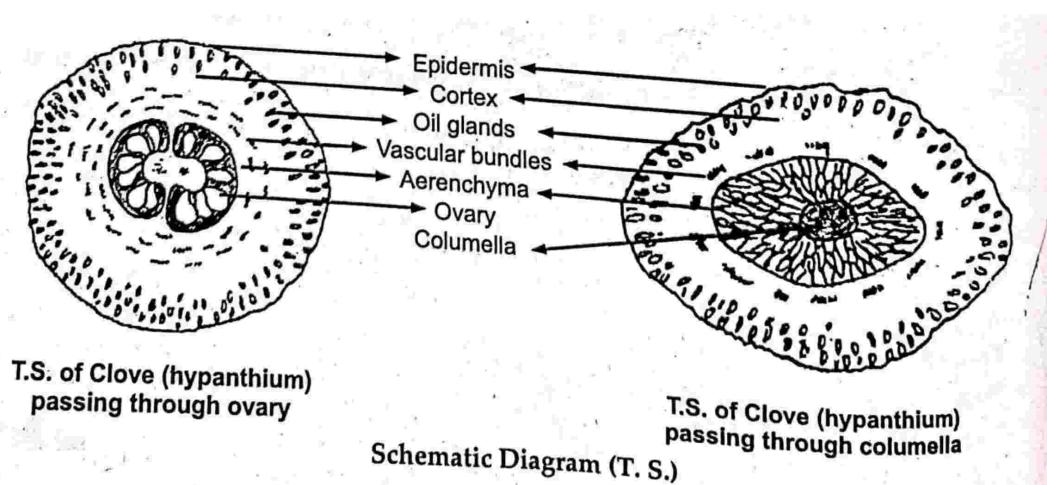




# Bright medico

(Never stop learning)

## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:

SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ HCl(1:1)	Conc. RED/PINK	Lignified tissues: xylem(vascular bundle))
2	Sudan Red III	red	Oil globule in the cells of endosperm and cuticle
3	Strong KOH solution	Needle shaped potassium eugenate crystals	Eugenol of Volatile oil
4	Dil hydrochloric acid	Crystals soluble	Calcium Oxalate crystals
5	Sulphuric acid	Soluble , needles of calcium sulphate on standing	Calcium oxalate crystals

### MICROSCOPICAL CHARACTERISTICS OF POWDERED DRUG :

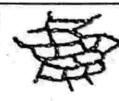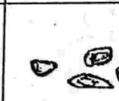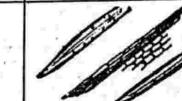Oil glands	Aerenchyma	Pollen grains	Fibres	Anther	Calcium oxalate crystals
					
Fragments of aerenchyma showing long oval schizolysigenous oil glands.	Portion of loose parenchyma with air spaces.	Diameter : 15 to 20 $\mu$ Small, biconvex, rounded or triangular in shape.	Sclerenchymatous fibres with parenchymatous cells	Fibrous layer of anther with reticulated cells	Sphaerophides No starch grains

### CHEMICAL TEST:

SL NO	TESTS	OBSERVATION	INTERFERENCE
1	Aqueous extract+ Lead acetate solution	White ppt	Tannins
2	Clove oil + alcohol+ ferric chloride 5% solution	Blue colouration	Eugenol
3	Aqueous extract + Ferric chloride solution 5%	Dark colour	tannins

**CHEMICAL CONSTITUENTS:** volatile oil: Eugenol, isoeugenol, methyl and dimethyl furfural, alpha and beta caryophylline , hydrolysable tannins.

**USES:** Carminative, aromatic, antiseptic, stimulant, flavouring agent, dental analgesic oil, in microscopic work , for isolation of eugenol.

# Bright medico


(Never stop learning)

## EXPERIMENT NUMBER -10

**AIM:** To study the macroscopic and microscopic study of the CINNAMOM


**Synonym:** Cinnamon bark, kalmi-dalchini

**Biological Source:** It Consists of dried inner bark of the shoots of coppiced trees of *Cinnamomum zelanicum*. It contains not less than 0.1% volatile oil. **Family:** Lauraceae.

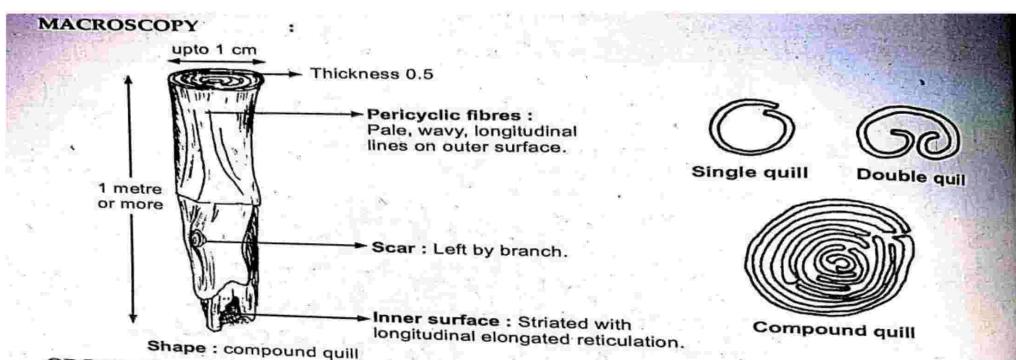

**ORGANOLEPTIC CHARACTERS:** Colour: outer surface dull yellowish brown, inner surface darker in colour, Odour: Fragrant, Taste: Warm, sweet.

**EXTRA FEATURES:** Bark is free of cork, single or double quills or single closely packed compound quill. **Fracture:** splintery.

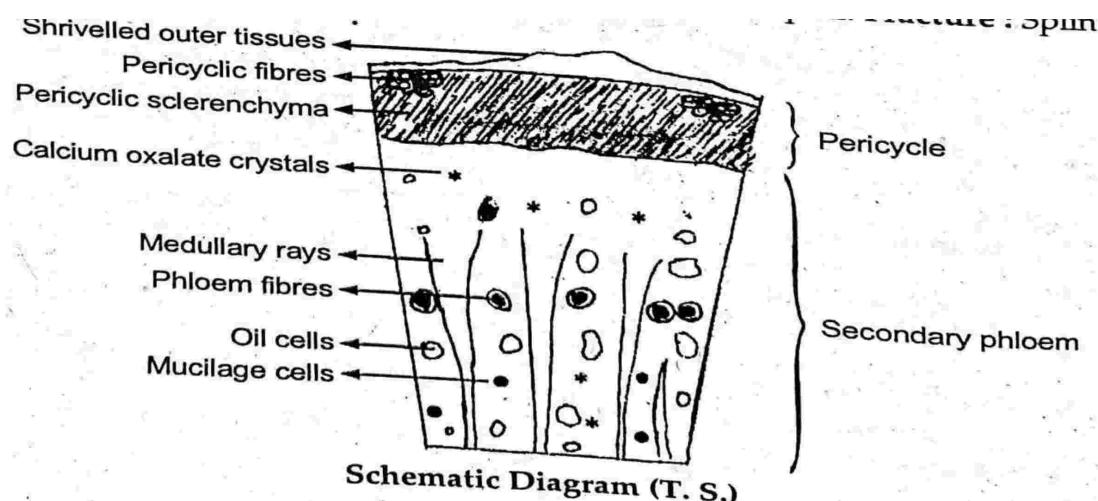


# Bright medico

(Never stop learning)


## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:

SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ Conc. HCl(1:1)	RED/PINK	Lignified tissues: xylem(vascular bundle))
2	Iodine	blue	strach
3	Ruthanium red	pink	Mucilage cells
4	Dil hydrochloric acid	Crystals soluble	Calcium Oxalate crystals
5	Acetic acid	insoluble	Calcium oxalate crystals
6	1% solution of osmic acid	Brown or pale	Volatile oil
7	Dil. Tincture alkana	Red on standing	Volatile oil


## MICROSCOPICAL CHARACTERISTICS OF POWDERED DRUG :

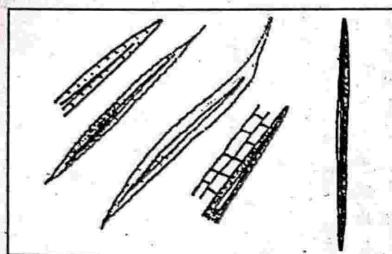


**Phloem Fibres :** Not more than 30  $\mu$  in diameter and 200 to 600  $\mu$  in length, stratified thick lignified wall and narrow lumen. The total length of fibres in cinnamon bark is 230 to 290 m per gram of air-dry bark.

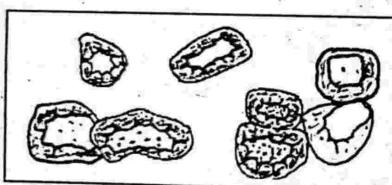

**Stone Cells :** U-shaped, lignified structures with one side flattened.

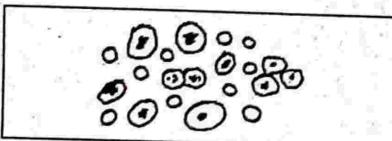

**Starch Grains :** Diameter below 10  $\mu$ ; present in parenchymatous cells of phloem and medullary rays.

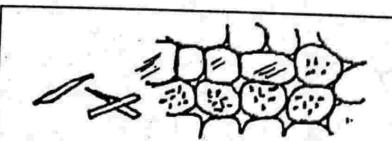

**Calcium Oxalate Crystals :** Small acicular raphides from parenchyma and medullary rays.

# Bright medico

(Never stop learning)

## CHEMICAL TEST:

SL NO	TESTS	OBSERVATION	INTERFERENCE
1	Volatile oil+ 5ml of alcohol+ one drop of ferric chloride	Green colour	Cinnamic aldehyde + eugenol present in volatile oil
2	Chloroform extract or volatile oil on slide+ drop of 10% aqueous phenyl hydrazine hydrochloride solution	Red shaped crystals	Cinnamic aldehyde
3	Aq. Solution + FeCl <sub>3</sub> solution	Dark colour	Tannin
4	Aq. Extract + lead acetate	White ppt	Tannin
5	AQ. Extract + potassium permanganate solution	Decolourisation	Tannin

**CHEMICAL CONSTITUENTS:** volatile oil 0.5-1%:Cinnamic aldehyde (55-65%) , Eugenol, tannins, starch , calcium oxalate.

**USES:** Carminative, aromatic, mild astringent, powerful germicide.


# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -11


**AIM:** To study the macroscopic and microscopic study of the  
**GINGER** Synonym: Sonth, zinger, Jamica ginger

**Biological Source:** It Consists of dried rhizome of *Zingiber officinale*, scrapped to remove the dark outer skin and dried in the sun. **Family:** Zingiberaceae.

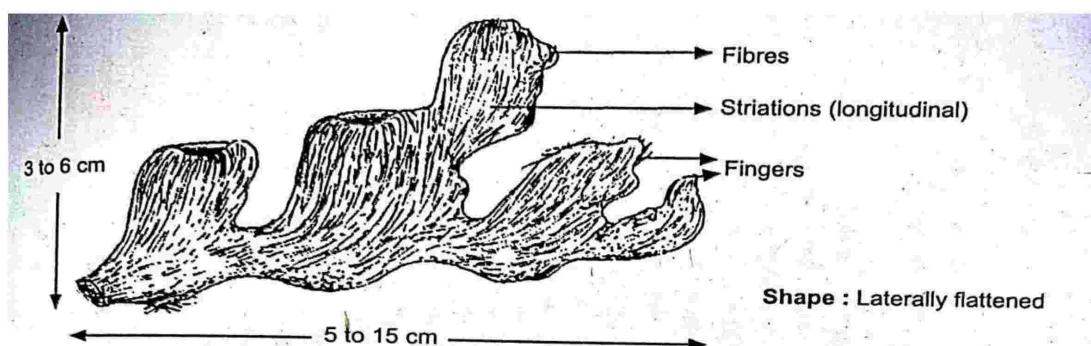

**ORGANOLEPTIC CHARACTERS:** Colour: externally buff colour, odour: agreeable and aromatic, Taste: agreeable, pungent, and characteristic.

**EXTRA FEATURE:** Sympodial branching, horizontal rhizome, transversely cut surface shows well marked endodermis and stele. **Fracture:** short and fibrous.

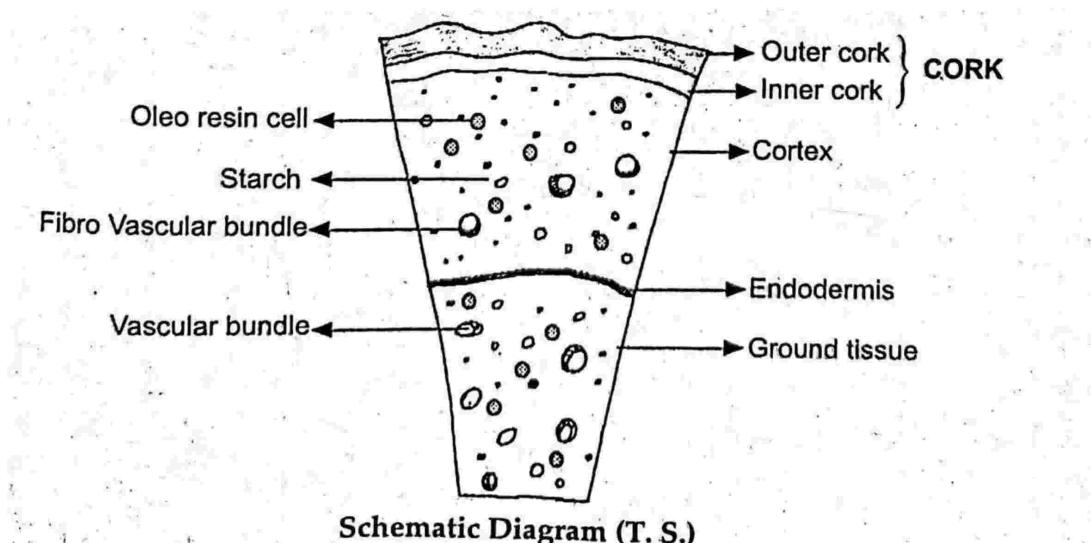


Schematic Diagram (T. S.)

# Bright medico

(Never stop learning)

## STRAINING/ DIAGNOSIS/MACRO-CHEMICAL TESTS:

SL NO	REAGENTS	OBSERVATION	CHARACTERISTICS
1.	Phloroglucinol+ HCl(1:1)	Conc. RED/PINK	Lignified tissues: xylem(vascular bundle))
2	Iodine	blue	starch

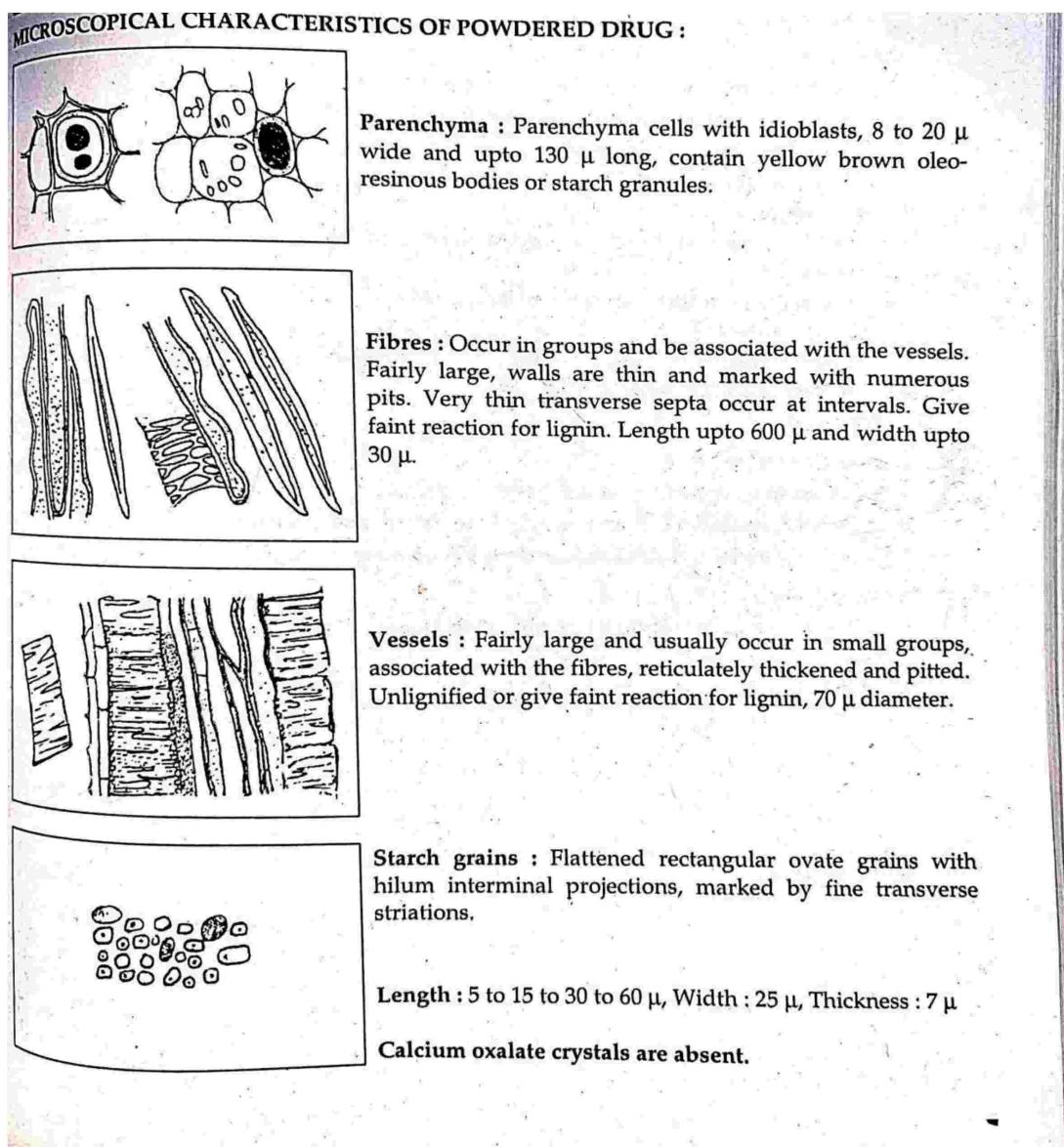

# Bright medico

(Never stop learning)

**CHEMICAL TEST:** Boil the drug with 5% potassium hydroxide or alkali – pungency of the ginger destroyed.

**CHEMICAL CONSTITUENTS:** volatile oil, Terpenes, cineol, citral, borneol, **pungent principle:** gingerol, shogaol, zingerone.

**USES:** Carminative, stimulant, flavouring agent

# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -12

**AIM:** To study the test of the AGAR:

**CHEMICAL REQUIREMENTS:** agar, ruthenium red, barium chloride, fehling solution, iodine solution.

**Test of Agar:**

Sl. NO	Experiments	Observation	Inference
1	Agar is boiled with water	Stiff jelly on cooling	Agar may be present
2	Agar solution is treated with ruthenium red	Pink color	Agar may be present
3	Hot agar solution is treated with barium chloride reagent	White PPT	Agar may be present
4	Hot agar solution is treated with fehling's solution and heated	Red PPT	Agar may be present
5	Agar solution is treated with iodine solution	Crimson to brown color	Agar may be present
6	Agar ash taken on a slide and add two drops HCl and observed under microscopes	Spongy spicules of diatoms are observed	Agar may be present

# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -13

**AIM:** To study the test of the acacia

**CHEMICAL REQUIREMENTS:** Acacia, Strong lead acetate, Sulphuric acid, Benzedrine, Iodine solution, HCl, Fehling solution

**Test for Acacia:**

SL NO	EXPERIMENTS	OBSERVATION	INFERENCE
1	5ml of a 2%W/V solution of acacia is treated with 1 ml of strong lead sub acetate solution	Flocculent white PPT is produced	Acacia may be present
2	Dissolve 0.25gm of acacia in 5 ml of water by shaking in cold. add 0.5 ml of H <sub>2</sub> SO <sub>4</sub> and 0.5 ml 1% solution of benzidine in alcohol, shake and allowed to stand	Unstable deep blue color is noticed (due to enzyme oxidase)	Acacia may be present
3	10ml of 2%W/V solution of acacia, add 0.2 ml of a 20%w/v solution of lead acetate	No ppt is produced	Acacia may be present
4	0.1 gm of acacia powder add 1 ml of N/50 iodine solution	The mixture does not acquire crimson color (distinction from agar and	Acacia may be present

# Bright medico

(Never stop learning)

		tragacanth)	
5	To 1 ml of acacia solution add 4 ml of water and dilute HCl acid and boiled for few minutes. Add Fehling's solution and heat	Red ppt (due to cuprous oxide)	Acacia may be present

# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -14

**AIM: To study the test of the TRAGACANTH**

**CHEMICAL REQUIREMENTS:** Tragacanth, Ferric chloride solution, Copper oxide, conc. Ammonium hydroxide, Lead acetate, Fehling solution, Caustic potash.

**Test for Tragacanth:**

Sl no	Experiment	Observation	Inferences
1	Drug is boiled with freshly prepared 10% aqueous ferric chloride solution	Deep yellow ppt	Tragacanth may be present
2	Dissolve tragacanth and precipitated copper oxide in conc. ammonium hydroxide	Stringy ppt	Tragacanth may be present
3	Tragacanth solution is treated with Fehling's solution and heated	Red ppt	Tragacanth may be present
4	Drug solution treated with lead acetate solution	White ppt	Tragacanth may be present
5	Tragacanth is treated with 5% aqueous caustic potash	Canary yellow color	Tragacanth may be present

# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -15

**AIM: To study the test of the GELATIN**

**Chemical requirements:**

- ✓ Gelatin
- ✓ Copper sulphate
- ✓ Conc. Nitric acid
- ✓ Mercuric nitrate
- ✓ Ninhydrin solution
- ✓ Nitrous acid

Sl. No	Test	Observation	Inference
1	Biuret test: To alkaline solution of a protein (2ml), a dilute solution of copper sulphate is added	A red or violet color is obtained with peptide containing at least two peptide linkage	
2	Xanthoproteic reaction: Protein + warmed with conc. Nitric acid	Give yellow fever	The color becomes orange color when solution made alkaline
3	Millon's reagent: (Mercuric nitrate in Nitric acid containing a trace of nitrous acid)	Give white PPT	PPT turn on red color on heating.
4	Aq. Solution of protein + alcoholic solution of Ninhydrin and heated	Red to violet color	

# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -16

**AIM: To study the test of the STARCH**

**Chemical requirements:**

- ✓ Starch
- ✓ Iodine solution

Sl. No	Test	Observation
1	Boil 1 gm of starch with 15 ml of water and cool	The translucent viscous jelly is produced
2	Adding iodine solution on the above jelly	Jelly turns deep blue in colour. The blue colour disappears on warming and reappear on cooling

# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -17

**AIM: To study the test of the HONEY**

**Chemical requirements:**

- ✓ Honey
- ✓ Fehling's solution A and B
- ✓ Molisch reagent ( $\alpha$ - Naphthol dissolved in Ethanol)
- ✓ Conc.  $H_2SO_4$

Sl. No	Test	Observation	Inference
1	Honey + Fehling solution A and B	Red PPT	Honey is present
2	Sugar solution + Molisch reagent +conc. $H_2SO_4$	Purple ring form	Honey is present

# Bright medico

(Never stop learning)

## EXPERIMENT NUMBER -18

**AIM: To study the test of the CASTOR**

**OIL Chemical requirements:**

- ✓ Petroleum Ether

Sl. No	Test	Observation
1	Castor oil treated with the half of its volume Petroleum Ether (50-60)	Completely soluble in petroleum ether
2	Oil +equal volume of alcohol and cool to 0C for 3 hrs.	Castor oil is present